Thanks for Choosing Zerhunt Oral Irrigator

Traditional Clean: Traditional tools can't clean hidden

food debris and harmful bacteria.

Zerhunt Oral Irrigator: Cleaning in places where traditional floss and toothbrushes can't reach.

Gently remove tiny food particles often trapped in the small

periodontal pockets between teeth and gums.



Contact us:

Email:zerhunt@outlook.com Zerhunt E-Business Co..Ltd




Portable Oral Irrigator Instruction Manual

FC159

Contents

Contents	-01
Product Description and Accessories	02
Getting Started ·····	.03
Instructions for Use	04
Recommended Technique	05
Remove and Clean Reservoir	06
Battery Maintenance and Disposal	.07
Troubleshooting Guide ·	08

Product Description and Accessories

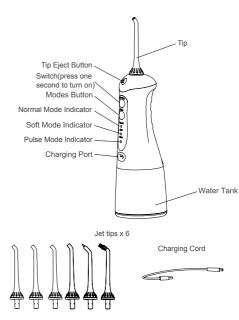

Charging

Turn off the unit and charge for 4 hours hefore first use Connect the charging cord with unit and plug the usb plug into your local standard outlet. The unit will not work when charging.



Insert tip until it clicks into place. and you will hear clicks.






Do not press tip eiect button while unit is running.



Lift the flip top on reservoir and fill warm water. Note: Please turn off unit while filling water.

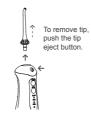

- . Press one second to turn on, don't worry about turning it on by touched accidentally when using or on your bag.
- . Comes with non-volatile memory function, turn power on, the mode is the same as the last time you used.
- . Mode switching is available when power off. You can choose the mode you want before turn on the unit.
- Press Mode button to switch between Normal. Soft and Pulse.
- . If one mode lasts for more than 120 seconds, the unit will turn off automatically.
- . Press Mode button for about 3s to indicate how much battery is left.
- . When charging, the 3 lights will flash in sequence. Till the 3 lights stay on, that means the unit has been fully charged.
- . In use, the 3 lights are flashing that means lacking of electricity, please recharge the unit.


Manually rotate the tip to change stream direction



Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth and turn the unit on.




Turn power on, close lips enough to prevent splashing. Keept unit upright during use to unsure continuous water flow. Aim water at the gumilne at a 90 degree angle. Follow the gumilne and pause briefly between teeth.

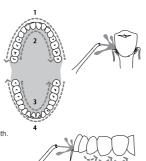

After Use

Turn power off. Empty liquid left in the water tank.



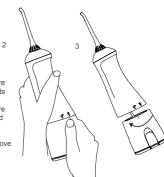

As shown in figure 2. One hand holds body of unit. As shown in figure 3. The other hand seize the water tank and turn clockwise to remove water tank

Battery Maintenance and Disposal

Battery Maintenance

- Please fully charge the unit before first use, it will take about
 hours.
- 2. The unit will not work when in charging.
- 3. Please fully charge the unit before storing it for a long time.
- To maximize the lifespan of the battery, please recharge the unit, when in low battery.

Battery Disposal

The unit contains a lithium battery, at the end of battery life, please put the unit at recycle center per your local or state requirements. The built-in battery can not be replaced.

Troubleshooting Guide

Problem	Cause	Solution
water leaks between tips and white body	improper installation of tip or gasket in tip base broken	make sure tip installed properly
water tank leaks	improper installation water tank broken O ring seal on flip damaged	Install water tank properly per instruction Contact Seller
inadequate pressure or can not turn on	low power	recharge the unit

If still can not solve the problems. Please contact seller to get after-sale support.

Specifications

◆ Power: 5w ◆ Adaptor: DC5V.1A

◆ Reservoir: 10.25oz / 300ml

◆ Pressure Range: 30~100psi◆ Pulse Rate: 1400~1800 times/min

♦ Mode: Normal, Soft, Pulse

◆ Battery: 1400mAh lithium battery